

medDream

PACS

**STORE, ARCHIVE,
MANAGE, VIEW
MEDICAL IMAGES**



VENDOR NEUTRAL
ARCHIVE



MODALITY
WORKLIST



CUSTOMISABLE
AUTOROUTING




HIS / RIS / EMR

HIS / RIS / EMR
INTEGRATION

 **Softneta**

Workflow with MedDream PACS

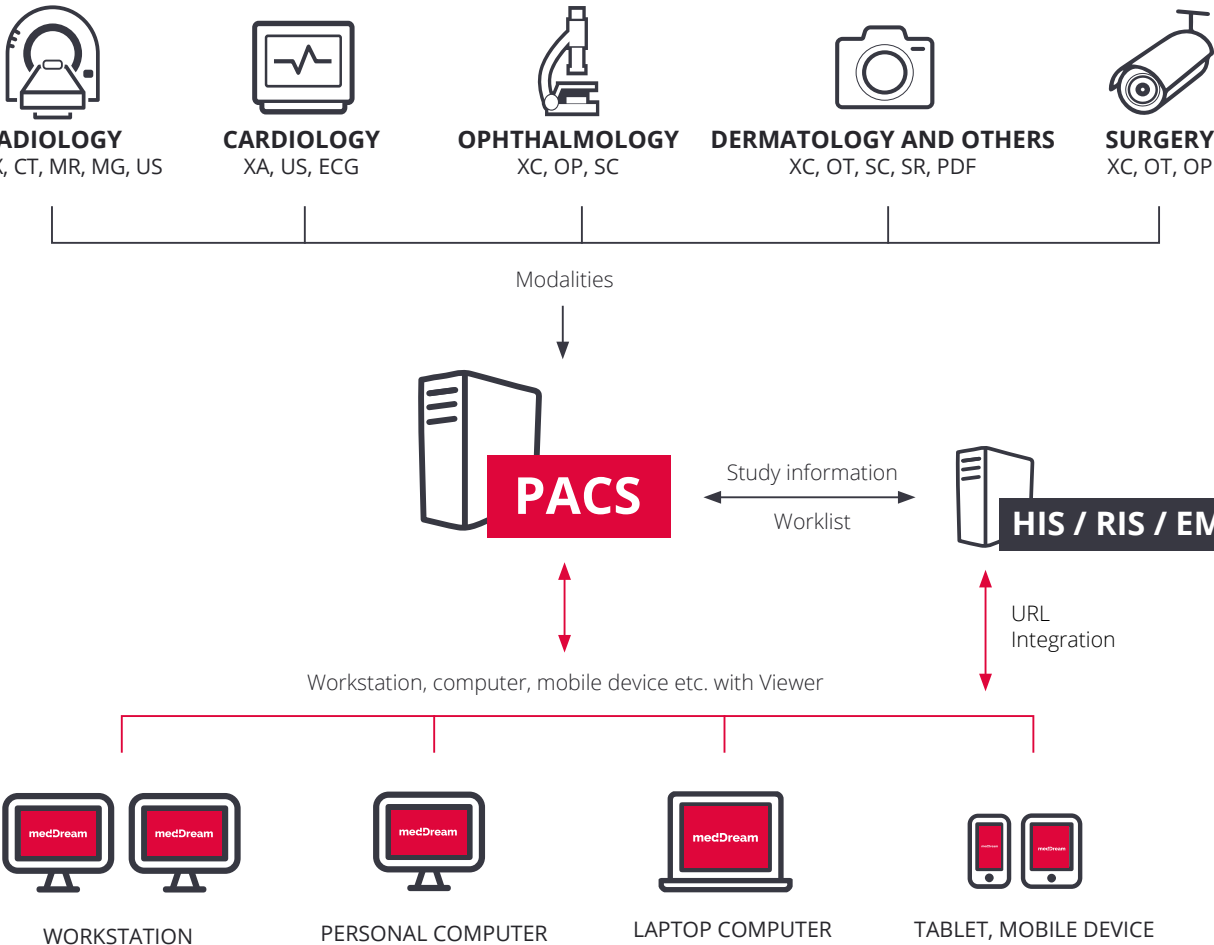

To	1d	3d	7d	1m	3m	Any	All	CR	CT	DX	MR	+	+	+	+
Patient ID	Name	Accessionnum	Modality	Description	Date Time	Received On	Source All								
A 1st0104	Anonymous	20160317145141	CT		2016-03-17 14:51:41	2016-03-17 14:23:06	MEDREAM								
A 1st01	Anonymous	20160318153612	CT		2016-03-18 15:36:12	2016-03-17 08:36:01	Upload default								
A 00-028	Anonymous	20160204182024	US		2016-02-04 18:20:24	2016-02-04 18:29:27	Upload root								
A 00-027	Anonymous	20160204142349	US		2016-02-04 14:23:49	2016-02-04 18:27:47	Upload root								
A 1st01	Anonymous	20160204101334	US		2016-02-04 10:13:34	2016-03-01 15:50:00	SENDTOPACS								
A 15-4353	Anonymous	2016020109305408	OT		2016-02-01 09:30:54	2016-02-04 11:43:28	Upload root								
A 123456	Anonymous	2016010815545289	OT		2016-01-08 15:54:52	2016-03-01 16:16:22	SENDTOPACS								
A 123456	Anonymous	2016010808000579	OT		2016-01-08 08:00:29	2016-03-01 16:18:13	SENDTOPACS								
A 1st01	Anonymous	5040	US	SendTOPACS DEMO	2015-12-09 13:42:48	2016-02-10 11:57:12	MEDREAM								
A CR15-180613	Anonymous	247-250	CT	KUNAS	2015-12-08 11:19:48	2016-01-26 11:51:51	MEDREAM								
A 1615-6648	Anonymous	182-185	CT	KUNAS	2015-12-07 07:31:53	2016-01-29 11:38:48	MEDREAM								
A 15-11528	Anonymous	201503201124348	OT		2015-03-23 11:24:33	2015-06-17 18:12:49	SENDTOPACS								
A 0	Anonymous	2014112712356688	OT		2016-11-27 12:56:44	2016-02-04 11:44:02	Upload root								
[14-198[SENDTOPACS-0411181225]	Anonymous	198	MR	GALVA+C	2014-06-19 14:32:38	2015-06-17 18:12:34	SENDTOPACS								
A [14-16	Anonymous		MR	GALVA+KALAS	2014-08-28 10:55:23	2015-06-17 18:11:54	MEDREAM								
A 4581118044	Anonymous		US		2014-02-12 11:38:37	2016-03-01 14:55:23	Upload default								
A 40315430521	Anonymous		MR	ABDOMEN+PELVIS	2014-01-29 10:09:31	2015-06-17 18:12:04	SENDTOPACS								
A 48411620312	Anonymous		MR	HEAD+ROUTINE	2014-04-28 08:52:18	2015-06-17 18:13:20	SENDTOPACS								
A 14-205-815	Anonymous	977-878	CT	STUBURAS 215	2014-01-27 07:52:11	2015-06-17 18:12:42	MEDREAM								
A 45403070645 - 4392	Anonymous		CT	Spine/Spine Routine (Subst)	2013-10-21 19:11:20	2015-06-17 18:13:27	MEDREAM								

1 to 20 of 89

Previous
1
2
3
4
5
Next

Copyright 2019 © Suttera MedStream 4.07.3 (Build 2019)

MedDream PACS + WEB DICOM Viewer features

Store, archive, manage medical images:

- C-STORE, C-FIND, C-GET and C-MOVE support;
- Implicit Little Endian, Explicit Little Endian;
- JPEG, JPEG-LS, JPEG 2000, RLE, MPEG-2;
- MPEG-4 Transfer Syntaxes supported;
- DICOM ECG file support;
- DICOM MPEG-2 and MPEG-4 Video playback support;
- Study forwarding.

View medical images:

- Regular tools: zoom, pan, windowing, magnifier, etc.;
- Cine playback of multi-frame sequences;
- Wide variety of measurements: distance, angle, density (HU), etc.;
- Reference lines (Scout Lines);
- Comfortable bar of series preview with thumbnails.

Integration:

- Integration to HIS / RIS / EMR. Flexible and open integration interface;
- Multi-language support.



Optional modules of MedDream WEB DICOM Viewer



Report Module

Writing, storing and viewing reports, templates





ECG Module

Viewing and measuring



Video Module

MPEG-2, MPEG-4 viewing

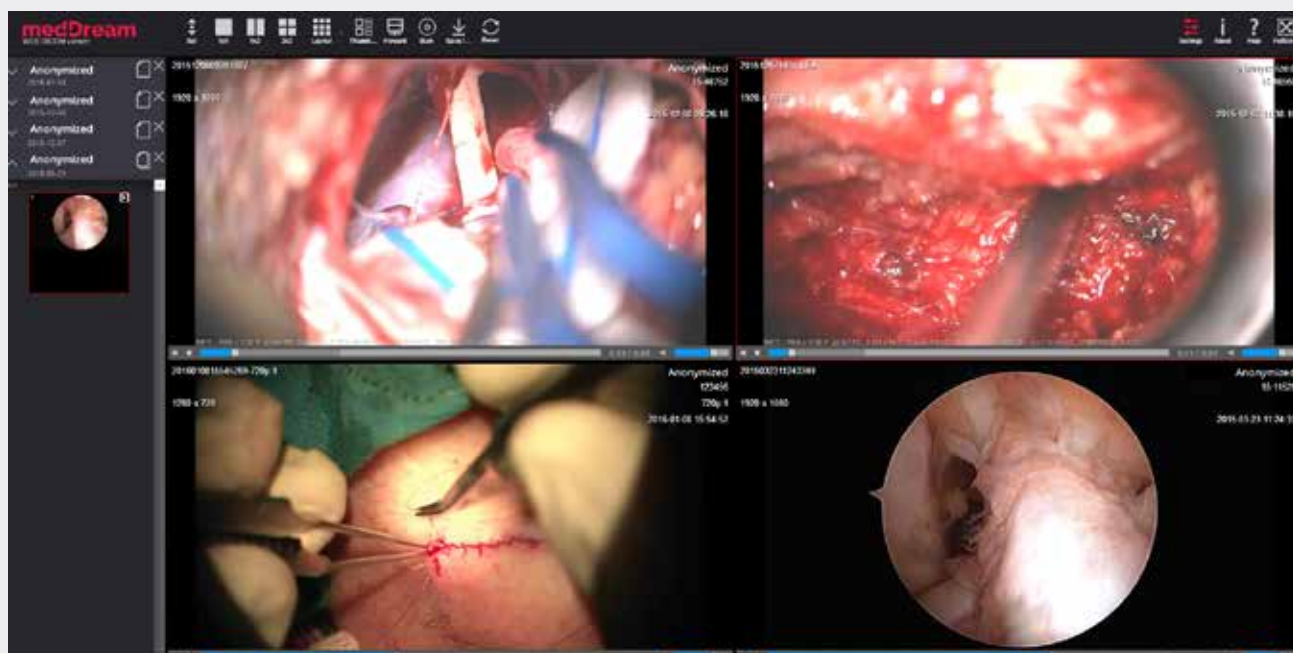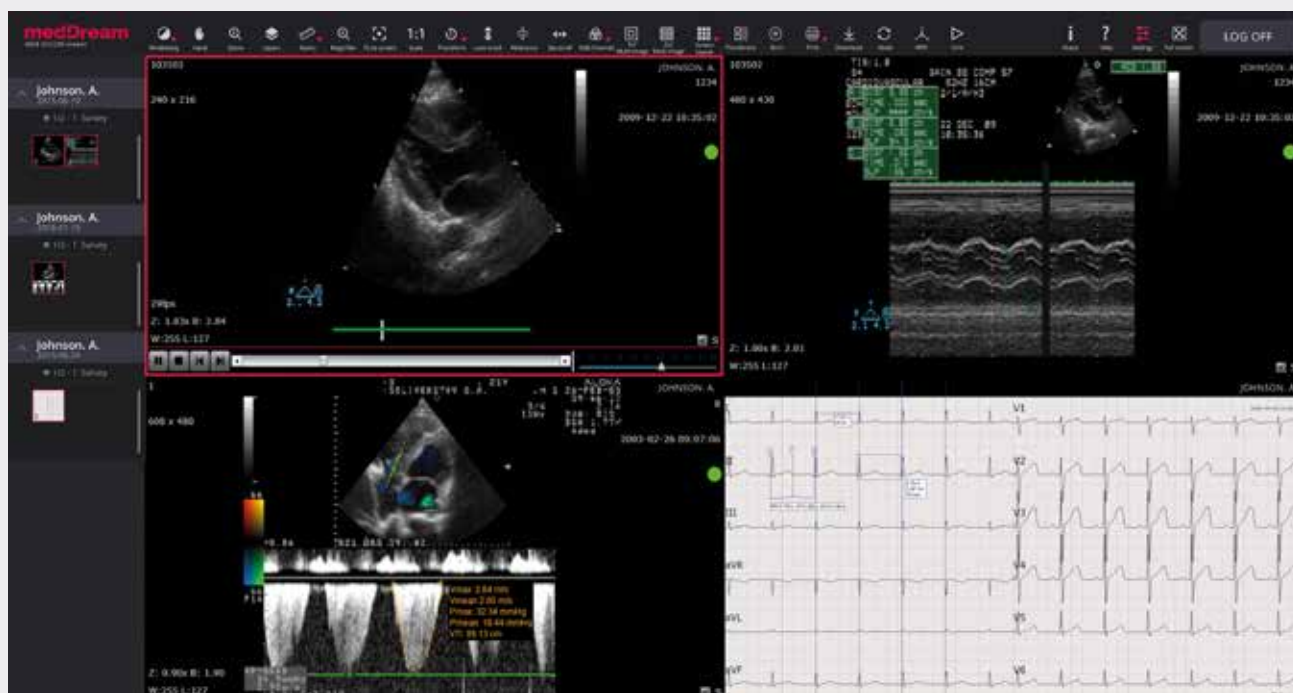

E-mail: info@lisit.jp
<http://www.lisit.jp>

東京都渋谷区円山町6-7
アムフラット1F

開発元 (Developed by)

